



SimplySnap **FEATURES**

Offset Electrodes



REDUCE ARTIFACTS

Offset connector minimizes artifact from patient & cable movements.



INTUITIVE APPLICATION

Clinicians can attach cables while electrodes are on patient without compromising gel well.



WET GEL

Wet gel provides quick reading & stabilization.

VersaTrode⁺ **FEATURES**

Snap 'n Tab Electrodes



ATTACHMENT FLEXIBILITY

Electrode can be used with either a snap or an alligator clip.



VERSATILE AND DURABLE

Suitable for all diagnostic procedures and stays on patient through multiple departments.

Compares to both the BlueSensor[®] and Skintact[®] Tab 'n Snap.

PACKAGING



RESEALABLE BULK POUCH

Ziplock pouch keeps electrodes fresh after opening.



CASE

FDA / GS1 Compliant Labels



1.800.669.6905 | HS.NisshaMedical.com

Patient Monitoring & Medical Surgical Consumables

MR Conditional



Not made with natural rubber latex

Wet Gel



VersaTrode⁺

Snap 'n Tab Electrodes

SimplySnap

Offset Electrodes

Description	Product #	Reference	Size	Snap	Per Pouch	Per Box	Icons	CROSS REF	
								SkinTact	Ambu
SimplySnap Offset Electrode	4010425C	OTL-5235-S	52 x 35	Stainless	60	1200		FS-VB01	Ambu R-00-S Ambu T-00-S Ambu SP-00-S
	32029077	OTL-5235-S-5	52 x 35	Stainless	5	1000			
SimplySnap Radiolucent Offset Electrode	4010452C	OTL-5235-SRT	52 x 35	Carbon	60	1200		N/A	N/A
	32029078	OTL-5235-SRT-5	52 x 35	Carbon	5	1000			
Offset Electrode No Snap	4010423C	OTL-5235	52 x 35	None	100	1500		FS-VBOO	Ambu SU-00-C
	32029079	OTL-5235-5	52 x 35	None	5	1000			
SimplySnap Pediatric Offset Electrode	4010453C	OTL-4227-S	42 x 27	Stainless	60	1200		F-VK01	Ambu P-00-S/50
	32029080	OTL-4227-S-5	42 x 27	Stainless	5	1000			
VersaTrode+ Snap 'n Tab	4010429C	TSL-4532	45 x 32	Carbon	60	1200		FS-WA10C/5	Ambu BS-Tab-2300/100
	32029076	TSL-4532-5	45 x 32	Carbon	5	1000			

ANATOMY



OFFSET SNAP

Offset snap allows clinician to place electrode on the patient first and then attach cable.




FOAM SUBSTRATE

Comfortable conforming foam provides strength without stretching.

SimplySnap

Offset Electrodes

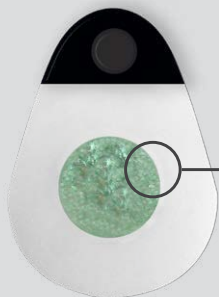


WET GEL

Wet gel provides quick reading and stabilization for faster traces.

SNAP OR TAB CONNECTION

Our most versatile electrode, allows connection with either a snap or a tab cable. Provides the most flexibility for use in all departments.

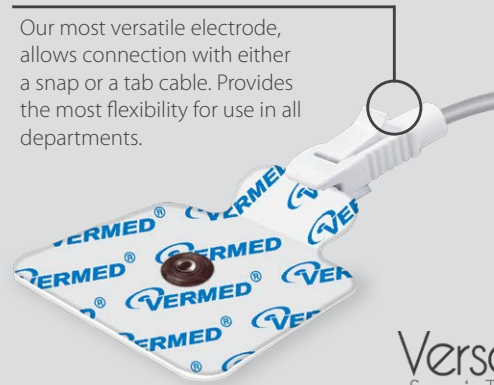

VersaTrode⁺

Snap 'n Tab Electrodes



ADHESIVE

Made from medical grade acrylic - ensuring that your electrodes stay in place.